## MT 200

IB MT 200 V13-2 0811

**Description of this Thermometer**

- ON/OFF Button
- Display
- Battery Compartment Cover
- Measuring Sensor

**Important Safety Instructions**

- The instrument may be used only for measuring body temperature!
- The minimum measurement time until the beep is heard must be maintained without exception!
- Ensure that children do not use the instrument unsupervised; some parts are small enough to be swallowed.
- Protect the instrument from impact and dropping!
- Avoid bending the tip more than 45°!
- Avoid ambient temperatures above 60 °C. NEVER boil the instrument!
- Use only the commercial disinfectants listed in the section «Cleaning and Disinfecting» to clean the instrument. The instrument must be intact when immersed in liquid disinfectant.
- We recommend this instrument is tested for accuracy every two years or after mechanical impact (e.g. being dropped). Please contact Microlife-Service to arrange the test.

 Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

-  Read the instructions carefully before using this device.

 Type BF applied part

**Turning on the Thermometer**

To turn on the thermometer, press the ON/OFF button 1; a short beep signals «thermometer ON». A display test is performed. All segments should be displayed.

The last measurement reading will be shown on the display 2 automatically for 2 seconds with the **«M»** icon.

Then at an ambient temperature of less than 32 °C, an **«L»** and a flashing **«°C»** appear at the display field 2. The thermometer is now ready for use.

**Function Test**

Correct functioning of the thermometer is tested automatically each time it is turned on. If a malfunction is detected (measurement inaccuracy), this is indicated by **«ERR»** on the display, and a measurement becomes impossible. In this case, the thermometer must be replaced.

**Using the Thermometer**

**IMPORTANT:** Place the thermometer at the measuring site before switching it on in order to achieve the most accurate result.

Choose the preferred measuring method. When taking a measurement, the current temperature is continuously displayed and the **«°C»** symbol flashes. If the beep is heard 10 times and the **«°C»** is no longer flashing, the predictive end-temperature has been determined and the thermometer can be read now.

10 short beeps will sound when the temperature is higher than 37.5 °C in order to alert the patient that he/she may have fever.

To prolong the battery life, turn off the thermometer by briefly pressing the ON/OFF button 1. Otherwise the thermometer will automatically turn off after about 10 minutes.

**Measuring Methods**

► **In the mouth (oral)**

Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor 4 must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/ exhaled air.

**Approx. measuring time: 10 seconds!**

► **In the anus (rectal)**

This is the most reliable measuring method, and is especially suitable for infants and small children. Carefully insert the measuring sensor 4 of the thermometer 2 to 3 cm into the anal aperture.

**Approx. measuring time: 10 seconds!**

► **In the armpit (axillary)**

To receive more reliable results we recommend to measure temperature orally or rectally.

**Cleaning and Disinfecting**

Name: Isopropyl alcohol 70 %; immerse: max. 24 hours.

**Battery Replacement**

When the **«▼»** symbol (upside-down triangle) appears at the display field, the battery is flat and needs replacing. To replace the battery remove the battery compartment cover 3 from the thermometer. Insert the new battery with the + at the top. Make sure you have a battery of the same type to hand. Batteries can be purchased at any electrical store.

**Technical Specifications**

<b>Type:</b>	Predictive maximum thermometer
<b>Measurement range:</b>	32,0 <span> </span> °C to 42,9 <span> </span> °C <p>Temp. &lt; 32,0<span> </span>°C: display <b>«L»</b> for low (too low) Temp. &gt; 42,9<span> </span>°C: display <b>«H»</b> for high (too high)</p>
<b>Measurement accuracy:</b>	± 0,1 <span> </span> °C between 34 <span> </span> °C and 42 <span> </span> °C
<b>Operating humidity:</b>	10 <span> </span> % to 40 <span> </span> %, 10-95 <span> </span> % relative maximum humidity
<b>Storage temperature:</b>	-25 <span> </span> °C to +60 <span> </span> °C; 10-95 <span> </span> % relative maximum humidity
<b>Battery:</b>	1.5/1.55 V; SR 41
<b>Reference to standards:</b>	EN 12470-3, clinical thermometers; ASTM E1112; IEC 60601-1; IEC 60601-1-2 (EMC)

This device complies with the requirements of the Medical Device Directive 93/42/EEC.

Technical alterations reserved.

**Guarantee**

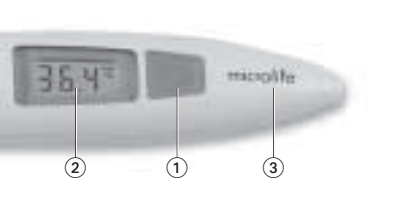
We grant you a lifetime guarantee after the date of purchase. Any damage caused by improper handling shall not be covered by the guarantee. Batteries and packaging are also excluded from the guarantee. All other damage claims excluded. A guarantee claim must be submitted with the purchase receipt. Please pack your defective instrument well and send with sufficient postage to the Microlife distributor.

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