# microlife®

## BP B2 Basic Blood Pressure Monitor

EN → 1	CZ → 35	HU → 67
RU → 9	SK → 43	HR → 75
BG → 18	SL → 51	PL → 83
RO → 27	SR → 59	

IB BP B2 Basic E-V11 0524  
Revision Date: 2024-01-02

Microlife Corporation  
9F, 431, RuiGuang Road, NeiHu  
Taipei 11492, Taiwan, China  
www.microlife.com

 Microlife UAB  
 P. Lukšio g. 32,  
08222 Vilnius  
Lithuania

CE0044

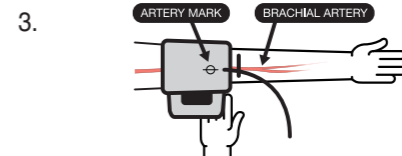
### Preparation

Sit on a back-supported chair and keep your legs uncrossed. / Usiądź na krześle wspieranym plecami i nie rozstawiaj nóg. / Сядьте, не скрещивая ноги, на стул со спинкой. / Сяднете на стол с облегалка и не кръстосвайте краката си. / Așezați-vă într-un scaun cu spatar și nu încrucișați picioarele / Sedněte si na židli s opěrkou a nepřekřičujte nohy. / Sadnite si na stoličku s opierkou a neprekrížujte si nohy. / Usedite se na stol s podprtim hrbtom in nogami na tleh (ne prekrížajte jih). / Sedite na stolicu koja podupire leđa i nemojte prekrštati noge. / Ülön egy háttámlás székre és lábait ne rakja keresztbe! / Sjednite na stolicu koja podupire leđa i neka Vam noge ne budu prekrizene.

Avoid thick or close-fitting garments on the upper arm. / Unikaj grubej lub dopasowanej odzieży na ramieniu. / Освободите плечо от плотной или плотно облегающей одежды. / Избягвайте дебели или плътно прилепващи дрехи върху горната част на ръката. / Evitați îmbrăcăminte groasă sau stâmta pe braț / Vyhněte se těsnému nebo přílišnému oděvu na horní části paže. / Usedite se na stol s podprtim hrbtom in nogami na tleh (ne prekrížajte jih). / Sedite na stolicu koja podupire leđa i nemojte prekrštati noge. / Ülön egy háttámlás székre és lábait ne rakja keresztbe! / Sjednite na stolicu koja podupire leđa i neka Vam noge ne budu prekrizene.

Place the artery-mark on the cuff over your artery. / Umieść mankietę na tętnicy nad arterią ręki w siłbie łokcia. / Поместите манжету так, чтобы значок артерии на манжете находился над артерией руки в си́бле локтя. / Поставете маркера на маншета върху артерията. / Plasați semnul artera al manșetei peste artera dumneavoastră / Umistěte manžetu se značkou tepny nad Vaší tepnu. / Umiestnite manžetu so značkou tepny nad Vašu tepnu. / Namestite znak za arterijo, ki je na manšeti, na vašo arterijo. / Postavite oznaku za arteriju na manžetni preko Vaše arterije. / Helyezze a mandzsetta artériajelzőjét az artériája fölé! / Postavite oznaku arterije na manžeti preko Vaše arterije.

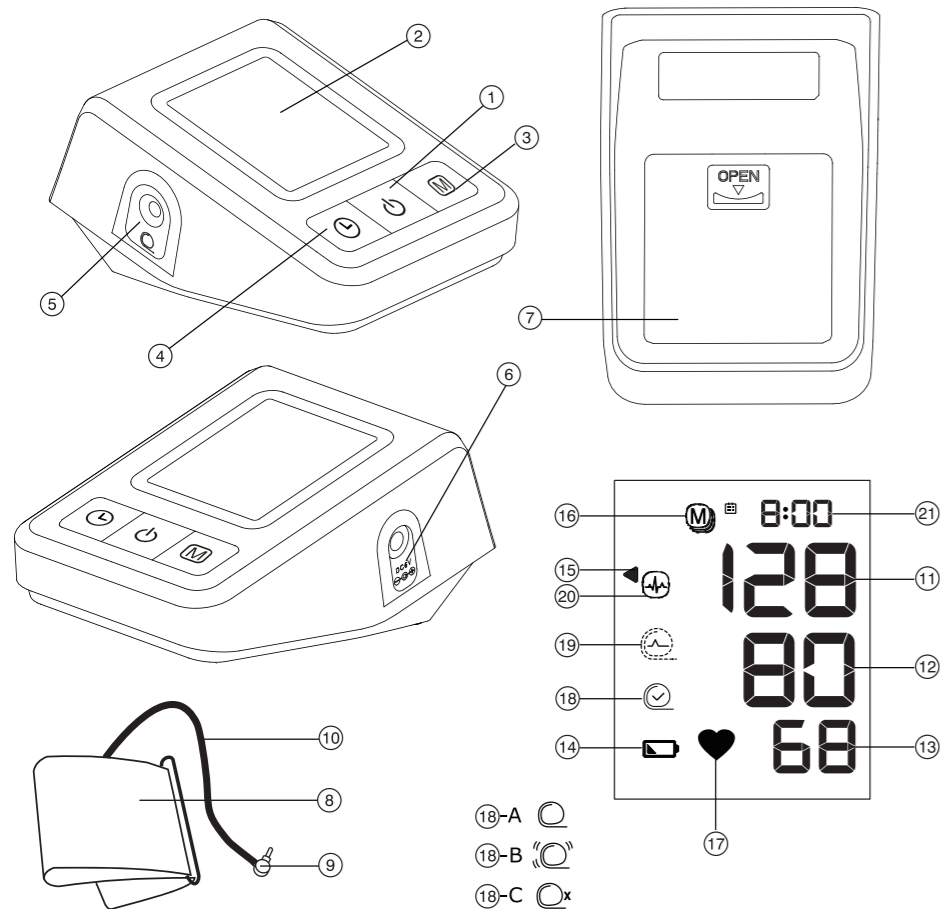
Fit the cuff closely, but not too tight. / Zacišnj mankiet dokladnie, jednak niezbyt silnie. / Поместите манжету так, чтобы она плотно, но не туго прилегала. / Стегнете добре маншета, но не твърде много. / Fixați manseta ferm dar nu prea stransă / Nasadte manžetu těsně, ale ne příliš těsně. / Nasadte manžetu tesne, ale nie príliš tesne. / Manšeto dobro namestite, vendar ne pretesno. / Dobro zategnite manžetnu, ali ne previše stegnuto. / Rögzítse a mandzsettát szorosan, de ne feszítse meg túl erősen! / Postavite manžetu blizu, ali ne pretjesno.

Position the cuff 1-2 cm above your elbow. / Ustaw mankiet na 1-2 cm powyżej łokcia. / Расположите манжету на 1-2 см выше локтя. / Поставете маншета 1-2 см над лакътя. / Plasați manșeta cu 1-2 cm deasupra articulației / Manžetu umistěte 1-2 cm nad loket. / Manžetu umiestnite 1-2 cm nad lakef. / Namestite manšeto 1-2 cm nad komolcem. / Postavite manžetnu 1-2 cm iznad Vašeg lakta. / A mandzsetta 1-2 centiméterrel legyen a könyökétől fejjebb! / Postavite manžetu 1-2 cm iznad Vašeg lakta.

Keep your arm still and do not speak during the measurement. / Trzymaj rękę nieruchomo i nie mów podczas pomiaru. / Держите руку неподвижно и не разговаривайте во время измерения. / Дръжте ръката си неподвижна и не говорете по време на измерването. / Mějte ruku uvolněnou a během měření nemluvte. / Maite ruku uvolnenú a počas merania nerozprávajte. / Držite ruku uvolnenu a počas merania nerozprávajte. / Med merjenjem ne premikajte roke in ne govorite. / Držite ruku mirno i nemojte pričati tokom merjenja. / Ne mozgassa a karját és ne beszéljen mérés közben! / Držite ruku mirno i nemojte pričati tijekom mjerenja.

1.

Avoid eating, bathing, smoking or caffeine (approx. 30 min).  
 Unikaj jedzenia, kąpieli, palenia lub kofeiny (około 30 minut).  
 Избегайте приема пищи, купания, курения или приема кофеиносодержащих продуктов (прибл. за 30 минут до измерения).  
 Избягайте ядене, къпане, пушене или кофеин (около 30 минути).  
 Evitajte sa mâncați, sa faceți baie, sa fumați, sau sa beți cafea (aproximativ 30 min).  
 Vyhněte se jídlu, koupání, kouření nebo kofeinu (přibližně 30 minut).  
 Vyhněte sa jedlu, kúpaniu, fajčeniu alebo kofeínu (približne 30 minút).  
 Izbegavajte se hrani, kupanju in kajenju (približno 30 minut).  
 Izbegavajte hranu, kupanje, pušenje i kofein (približno 30 minuta).  
 Kerülje az evést, fürdést, dohányzást vagy koffeinfogyasztást (a megelőző kb. 30 percben)!  
 Izbjegavajte jesti, kupati se, pušiti i kofein (oko 30 min.)

2.

Avoid activity and relax for 5-10 min.  
 Unikaj aktywności i zrelaksuj się przez 5-10 minut.  
 В течение 5-10 мин. оставайтесь в расслабленном состоянии и избегайте физической нагрузки.  
 Избягайте активности и се отпуснете в продължение на 5-10 минути.  
 Evitajte orice activitate si relaxati-vă 5-10 min.  
 Vyhněte se fyzické aktivitě a odpočívajte 5-10 minut.  
 Vyhnite sa fyzickej aktivite a odpočívajte 5-10 minút.  
 Izbegavajte se fizični aktivnosti in počivajte 5-10 minut.  
 Izbegavajte fizičku aktivnost i opustite se tokom 5-10 minuta.  
 Kerülje a fizikai tevékenységet, pihenjen 5-10 percig!  
 Izbjegavajte aktivnost i opustite se 5-10 min.

3.

Measure before medication intake.  
 Zmierz przed przyjęciem leku.  
 Измеряйте до приема лекарства.  
 Измерете преди приема на лекарството ви.  
 Măsurati înainte de a lua medicația.  
 Měření proveďte před přijetím léku.  
 Meranie uskutočnite pred užitím liekov.  
 Merjenje izvedite pred zaužitjem zdravil.  
 Obavite merenje pre uzimanja leka.  
 Mérjen a gyógyszer bevétele előtt!  
 Izvedite mjerjenje prije nego što uzmete lijek.

Name of Purchaser / Ф.И.О. покупателя /  
 Име на купувача / Numele cumpărătorului / Imię  
 i nazwisko nabywcy / Jméno kupujícího / Meno  
 zákazníka / Ime in priimek kupca / Ime i prezime  
 kupca / Vásártó neve / Ime i prezime kupca

Serial Number / Серийный номер / Серийен  
 номер / Număr de serie / Numer serijny /  
 Výrobní číslo / Výrobné číslo / Serijska številka /  
 Serijski broj / Sorozatszám / Serijski broj

Date of Purchase / Дата покупки / Дата на  
 закупуване / Data cumpărării / Data zakupu /  
 Datum nákupu / Dátum kúpy / Datum nakupa /  
 Datum kupovine / Vásárlás dátuma / Datum  
 kupovine

Specialist Dealer / Специализированный дилер /  
 Специалист дистрибутор / Distribuitor de spe-  
 cialitate / Przedstawiciel / Specializovaný dealer /  
 Špecializovaný predajca / Specializirani trgovec /  
 Ovlašćeni diler / Forgalmazó / Ovlašteni prodavač

- ① ON/OFF button
- ② Display
- ③ M-button (memory)
- ④ Time button
- ⑤ Cuff socket
- ⑥ Mains Adapter Socket
- ⑦ Battery compartment
- ⑧ Cuff
- ⑨ Cuff connector
- ⑩ Cuff tube

## Display

- ⑪ Systolic value
- ⑫ Diastolic value
- ⑬ Pulse rate
- ⑭ Battery display
- ⑮ Traffic light indicator
- ⑯ Stored value
- ⑰ Pulse indicator
- ⑱ Cuff fit check
  - A: Suboptimal cuff fit
  - B: Arm movement indicator «ERR 2»
  - C: Cuff pressure check «ERR 3»
- ⑲ Cuff signal indicator «ERR 1»
- ⑳ Irregular heartbeat (IHB) symbol
- ㉑ Date/Time

Read the important information in these instructions for use before using this device. Follow the instructions for use for your safety and keep it for future reference.

Type BF applied part

Keep dry

Manufacturer

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

Authorized representative in the European Community

Catalogue number

Serial number  
(YYYY-MM-DD-SSSSS;  
year-month-day-serial number)

Caution

Humidity limitation for operating **and** storage

Temperature limitation for operating **or** storage

Medical device

Keep away from children of age 0 - 3

CE Marking of Conformity

## Intended use:

This oscillometric blood pressure monitor is intended for measuring non-invasive blood pressure in people aged 12 years or older.

It is clinically validated in patients with hypertension, hypotension, diabetes, pregnancy, pre-eclampsia, atherosclerosis, end-stage renal disease, obesity and the elderly.

Dear Customer,

This device was developed in collaboration with physicians and clinical tests carried out prove its measurement accuracy to be of a very high standard.\*

If you have any questions, problems or want to order spare parts please contact your local Microlife-Customer Service. Your dealer

or pharmacy will be able to give you the address of the Microlife dealer in your country. Alternatively, visit the internet at [www.microlife.com](http://www.microlife.com) where you will find a wealth of invaluable information on our products.

Stay healthy – Microlife Corporation!

*\* This device uses the same measuring technology as the award winning «BP 3BTO-A» model tested according to the British and Irish Hypertension Society (BIHS) protocol.*

## Table of contents

### 1. Using the device for the first time

Inserting the batteries  
Setting the date and time  
Selecting the correct cuff

### 2. Checklist for taking a reliable measurement

### 3. Taking a blood pressure measurement

Manual inflation  
How not to store a reading  
How do I evaluate my blood pressure  
Appearance of the irregular heartbeat (IHB) symbol

### 4. Data memory

Viewing the stored values  
Clearing all values

### 5. Battery indicator and battery change

Low battery  
Flat battery - replacement  
Which batteries and which procedure?  
Using rechargeable batteries

### 6. Using a mains adapter

### 7. Error Messages

### 8. Safety, care, accuracy test and disposal

Device care  
Cleaning the cuff  
Accuracy test  
Disposal

### 9. Guarantee

### 10. Technical Specifications

## 1. Using the device for the first time

### Inserting the batteries

After you have unpacked your device, first insert the batteries. The battery compartment ⑦ is on the bottom of the device. Insert the batteries (4 x 1.5 V, size AA), thereby observing the indicated polarity.

### Setting the date and time

1. After the new batteries are fitted, the year number flashes in the display. You can set the year by pressing the M-button ③. To confirm and then set the month, press the time button ④.
2. Press the M-button to set the month. Press the time button to confirm and then set the day.
3. Follow the instructions above to set the day, hour and minutes.
4. Once you have set the minutes and pressed the time button, the date and time are set and the time is displayed.
5. If you want to change the date and time, press and hold the time button for approx. 3 seconds until the year number starts to flash. Now you can enter the new values as described above.

### Selecting the correct cuff

Microlife offers different cuff sizes. Select the cuff size to match the circumference of your upper arms (measured by close fitting in the centre of the upper arm).

Cuff size	for circumference of upper arm
S	17 - 22 cm
M	22 - 32 cm
M - L	22 - 42 cm
L	32 - 42 cm
L - XL	32 - 52 cm

 Pre-shaped cuffs are optionally available.

 Only use Microlife cuffs.

- ▶ Contact your local Microlife Service if the enclosed cuff ⑧ does not fit.
- ▶ Connect the cuff to the device by inserting the cuff connector ⑨ into the cuff socket ⑤ as far as it will go.

## 2. Checklist for taking a reliable measurement

- ▶ Avoid activity, eating or smoking immediately before the measurement.
- ▶ Sit down on a back-supported chair and relax for 5 minutes. Keep your feet flat on the floor and do not cross your legs.
- ▶ **Always measure on the same arm** (normally left). It is recommended that doctors perform double arm measurements on a patient's first visit in order to determine which arm to measure in the future. The arm with the higher blood pressure should be measured.
- ▶ Remove close-fitting garments from the upper arm. To avoid constriction, shirt sleeves should not be rolled up - they do not interfere with the cuff if they are laid flat.
- ▶ Always ensure that the correct cuff size is used (marking on the cuff).
  - Fit the cuff closely, but not too tight.
  - Make sure that the cuff is positioned 1-2 cm above the elbow.
  - The **artery mark** on the cuff (ca.3 cm long bar) must lie over the artery which runs down the inner side of the arm.
  - Support your arm so it is relaxed.
  - Ensure that the cuff is at the same height as your heart.

## 3. Taking a blood pressure measurement

1. Press the ON/OFF button ① to start the measurement.
2. The cuff will now pump up automatically. Relax, do not move and do not tense your arm muscles until the measurement result is displayed. Breathe normally and do not talk.
3. The cuff fit check ⑱ on the display indicates that the cuff is perfectly placed. If the icon ⑱-A appears, the cuff is fitted suboptimally, but it is still ok to measure.
4. When the correct pressure is reached, the pumping stops and the pressure falls gradually. If the required pressure was not reached, the device will automatically pump some more air into the cuff.
5. During the measurement, the pulse indicator ⑲ flashes in the display.
6. The result, comprising the systolic ⑩ and the diastolic ⑪ blood pressure and the pulse rate ⑬ are displayed. Note also the explanations on further display symbols in this booklet.
7. When the device has finished measuring, remove the cuff.

8. Switch off the device. (The monitor does switch off automatically after approx. 1 min.).

- ☞ You can stop the measurement at any time by pressing the ON/OFF button or open the cuff (e.g. if you feel uneasy or an unpleasant pressure sensation).
- ☞ This monitor is specially tested for use in pregnancy and pre-eclampsia. When you detect unusual high readings in pregnancy, you should measure after a short while again (eg. 1 hour). If the reading is still too high, consult your doctor or gynecologist.

### Manual inflation

**In case of high systolic blood pressure (e.g. above 135 mmHg)**, it can be an advantage to set the pressure individually. Press the ON/OFF button after the monitor has been pumped up to a level of approx. 30 mmHg (shown on the display). Keep the button pressed until the pressure is about 40 mmHg above the expected systolic value – then release the button.

### How not to store a reading

As soon as the reading is displayed press and hold the ON/OFF button ① until «M» ⑯ is flashing. Confirm to delete the reading by pressing the time button ④.

☞ «CL» is displayed when the reading is deleted from the memory successfully.

### How do I evaluate my blood pressure

The triangle on the left-hand edge of the display ⑮ points at the range within which the measured blood pressure value lies. The value is either within the optimum (white), elevated (hatched gray) or high (black) range. The classification corresponds to the following ranges defined by international guidelines (ESH, ESC, JSH). Data in mmHg.

Range	Systolic	Diastolic	Recommendation
1. blood pressure too high	≥135	≥85	Seek medical advice
2. blood pressure elevated	130 - 134	80 - 84	Self-check
3. blood pressure normal	<130	<80	Self-check

The higher value is the one that determines the evaluation.  
Example: a blood pressure value of **140/80** mmHg or a value of **130/90** mmHg indicates «blood pressure too high».

### Appearance of the irregular heartbeat (IHB) symbol

This symbol  indicates that an irregular heartbeat was detected. In this case, the measured blood pressure may deviate from your actual blood pressure values. It is recommended to repeat the measurement.

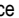
#### Information for the doctor in case of repeated appearance of the IHB symbol:

This device is an oscillometric blood pressure monitor that also measures the pulse during blood pressure measurement and indicates when the heart rate is irregular.

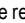
## 4. Data memory

This device automatically stores the last 30 measurement values.

### Viewing the stored values

Press the M-button  briefly, when the device is switched off. The display first shows «**M**» , and «**A**» which stands for the average of all stored values.

Pressing the M-button again displays the previous value. Pressing the M-button repeatedly enables you to move from one stored value to another.

-  Blood pressure readings with suboptimal cuff fit -A are not considered in the average value.
-  Pay attention that the maximum memory capacity of 30 memories is not exceeded. **When the 30 memory is full, the oldest value is automatically overwritten with the 31 value.** Values should be evaluated by a doctor before the memory capacity is reached – otherwise data will be lost.

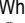
### Clearing all values

If you are sure that you want to permanently remove all stored values, hold down the M-button (the device must have been switched off beforehand) until «**CL ALL**» appears and then release the button. To permanently clear the memory, press the time button while «**CL ALL**» is flashing. **Individual values cannot be cleared.**

-  **Cancel deletion:** press ON/OFF button  while «**CL ALL**» is flashing.

## 5. Battery indicator and battery change

### Low battery

When the batteries are approximately  $\frac{3}{4}$  empty the battery symbol  will flash as soon as the device is switched on (partly filled battery displayed). Although the device will continue to measure reliably, you should obtain replacement batteries.

### Flat battery - replacement

When the batteries are flat, the battery symbol  will flash as soon as the device is switched on (flat battery displayed). You cannot take any further measurements and must replace the batteries.

- Open the battery compartment  at the back of the device.
- Replace the batteries – ensure correct polarity as shown by the symbols in the compartment.
- To set date and time, follow the procedure described in Section «1. Using the device for the first time».

-  The memory retains all values although date and time must be reset – the year number therefore flashes automatically after the batteries are replaced.

### Which batteries and which procedure?

-  Use 4 new, long-life 1.5 V, size AA alkaline batteries.
-  Do not use batteries beyond their date of expiry.
-  Remove batteries if the device is not going to be used for a prolonged period.

### Using rechargeable batteries

You can also operate this device using rechargeable batteries.

-  Only use «NiMH» type reusable batteries.
-  Batteries must be removed and recharged when the flat battery symbol appears. They should not remain inside the device as they may become damaged (total discharge as a result of low use of the device, even when switched off).
-  Always remove the rechargeable batteries if you do not intend to use the device for a week or more.
-  Batteries cannot be charged in the blood pressure monitor. Recharge batteries in an external charger and observe the information regarding charging, care and durability.

## 6. Using a mains adapter

You can operate this device using the Microlife mains adapter (DC 6V, 600 mA).

 Only use the Microlife mains adapter available as an original accessory appropriate for your supply voltage.

 Ensure that neither the mains adapter nor the cable are damaged.

1. Plug the adapter cable into the mains adapter socket  in the blood pressure monitor.

2. Plug the adapter plug into the wall socket.

When the mains adapter is connected, no battery current is consumed.

## 7. Error Messages

If an error occurs during the measurement, the measurement is interrupted and an error message, e.g. «ERR 3», is displayed.

Error	Description	Potential cause and remedy
«ERR 1» 	Signal too weak	The pulse signals on the cuff are too weak. Re-position the cuff and repeat the measurement.*
«ERR 2» 	Error signal	During the measurement, error signals were detected by the cuff, caused for instance by movement or muscle tension. Repeat the measurement, keeping your arm still.
«ERR 3» 	Abnormal cuff pressure	An adequate pressure cannot be generated in the cuff. A leak may have occurred. Check that the cuff is correctly connected and is not too loose. Replace the batteries if necessary. Repeat the measurement.
«ERR 5»	Abnormal result	The measuring signals are inaccurate and no result can therefore be displayed. Read through the checklist for taking a reliable measurement and then repeat the measurement.*
«HI»	Pulse or cuff pressure too high	The pressure in the cuff is too high (over 299 mmHg) OR the pulse is too high (over 200 beats per minute). Relax for 5 minutes and repeat the measurement.*

Error	Description	Potential cause and remedy
«LO»	Pulse too low	The pulse is too low (less than 40 beats per minute). Repeat the measurement.*

\* Please immediately consult your doctor, if this or any other problem occurs repeatedly.

## 8. Safety, care, accuracy test and disposal

### Safety and protection

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device may only be used for the purposes described in these instructions. The manufacturer cannot be held liable for damage caused by incorrect application.
- This device comprises sensitive components and must be treated with caution. Observe the storage and operating conditions described in the «Technical Specifications» section.
- The cuffs are sensitive and must be handled with care.
- Only pump up the cuff once fitted.
- Do not use this device if you think it is damaged or notice anything unusual.
- Never open this device.
- Read the additional safety information provided within the individual sections of this instruction manual.
- The measurement results given by this device is not a diagnosis. It is not replacing the need for the consultation of a physician, especially if not matching the patient's symptoms. Do not rely on the measurement result only, always consider other potentially occurring symptoms and the patient's feedback. Calling a doctor or an ambulance is advised if needed.

Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.

### Contra-indications

Do not use this device if the patient's condition meets the following contra-indications, to avoid inaccurate measurements or injuries.

- The device is not intended for measuring blood pressure in pediatric patients of age younger than 12 years old (children, infant, or neonates).
- Presence of significant cardiac arrhythmia during measurement may interfere with blood pressure measurement and affect the reliability of blood pressure readings. Consult with your doctor about whether the device is suitable for use in this case.
- The device measures blood pressure using a pressured cuff. If the measuring limb suffers from injuries (for example open wounds) or under conditions or treatments (for example intravenous drip) making it unsuitable for surface contact or pressurization, do not use the device, to avoid worsening of the injuries or conditions.
- Patient motions during measurement may interfere with the measurement process and influence results.
- Avoid taking measurements of patients with conditions, diseases, and susceptible to environment conditions that lead to uncontrollable motions (e.g. trembling or shivering) and inability to communicate clearly (for example children and unconscious patients).
- The device uses oscillometric method to determine blood pressure. The arm being measure should have normal perfusion. The device is not intended to be used on a limb with restricted or impaired blood circulation. If you suffer with perfusion or blood disorders, consult your doctor before using the device.
- Avoid taking measurement on the arm on the side of a mastectomy or lymph node clearance.
- Do not use this device in a moving vehicle (for example in a car or on an aircraft).

### **WARNING**

Indicates a potentially hazardous situation, which if not avoided, could result in death or serious injury.

- This device may only be used for the intended uses described in this Instructions for Use. The manufacturer cannot be held liable for damage caused by incorrect application.
- Do not change the patient medication and treatment based the result of one or multiple measurements. Treatment and medication changes should be prescribed only by a medical professional.
- Inspect the device, cuff, and other parts for damage. **DO NOT USE** the device, cuff or parts if they appear damaged or operating abnormally.

- Blood flow of the arm is temporarily interrupted during measurement. Extended interruption of blood flow reduces peripheral circulation and may cause tissue injury. Beware of signs (for example tissue discoloration) of impeded peripheral circulation if taking measurements continuously or for an extended period of time.
- Prolonged exposure to cuff pressure will reduce peripheral perfusion and may lead to injury. Avoid situations of extended cuff pressurization beyond normal measurements. In the case of abnormally long pressurization, abort the measurement or loose the cuff to depressurize the cuff.
- Do not use this device in oxygen rich environment or near flammable gas.
- The device is not water resistant or water proof. Do not spill or immerse the device in water or other liquids.
- Do not disassemble or attempt to service the device, accessory and parts, during use or in storage. Access to the device internal hardware and software is prohibited. Unauthorized access and servicing of the device, during use or in storage, may compromise the safety and performance of the device.
- Keep the device away from children and people incapable of operating the device. Beware of the risks of accidental ingestion of small parts and of strangulation with the cables and tubes of this device and accessories.

### **CAUTION**

Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury to the user or patient, or cause damage to the device or other property.

- The device is intended only for measuring blood pressure at upper arm. Do not measure other sites because the reading does not reflect your blood pressure accurately.
- After a measurement is completed, loosen the cuff and rest for > 5 minutes to restore limb perfusion, before taking another measurement.
- Do not use this device with other medical electrical (ME) equipment simultaneously. This may cause device malfunction or measurement inaccuracies.
- Do not use this device in proximity of high frequency (HF) surgical equipment, magnetic resonance imaging (MRI) equipment, and computerized tomography (CT) scanners. This may cause device malfunction and measurement inaccuracies.

- Use and store the device, cuff and parts in temperature and humidity conditions specified in the «Technical Specifications». Usage and storage of the device, cuff and parts in conditions outside ranges given in the «Technical Specifications» may result in device malfunction and the safety of usage.
- Protect the device and accessories from the following to avoid damaging the device:
  - water, other liquids, and moisture
  - extreme temperatures
  - impacts and vibrations
  - direct sunlight
  - contamination and dust
- Stop using this device and cuff and consult with your doctor if you experience skin irritation or discomfort.

### Electromagnetic Compatibility Information

This device is compliant with EN60601-1-2: 2015 Electromagnetic Disturbances standard.

This device is not certified to be used in vicinity of High Frequency (HF) medical equipment.

Do not use this device close to strong electromagnetic fields and portable radio frequency communication devices (for example microwave oven and mobile devices). Keep a minimum distance of 0.3 m from such devices when using this device.

### Device care

Clean the device only with a soft, dry cloth.

### Cleaning the cuff

The cuff delivered with this device is washable.

1. Remove the cuff connector (9) from the cuff tube (10) and carefully pull the bladder through the opening at the edge of the cuff cover.
2. Hand wash the cuff cover in soapsuds: not hotter than 30 °C.
3. Completely dry the cuff cover by linen drying.
4. Loop the cuff tube back through its opening and carefully place the bladder flat in the cuff cover.
5. Reattach the cuff connector on the cuff tube.

The bladder must lay straight in the cuff cover, not folded.

Do not use fabric softener.

**WARNING:** Do not wash the cuff in a washing machine or dishwasher!

**WARNING:** Do not dry the cuff cover in a tumble dryer!

**WARNING:** Under no circumstances may you wash the inner bladder!

### Accuracy test

We recommend this device is tested for accuracy every 2 years or after mechanical impact (e.g. being dropped). Please contact your local Microlife-Service to arrange the test (see foreword).

### Disposal

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

## 9. Guarantee

This device is covered by a **5 year guarantee** from the date of purchase. During this guarantee period, at our discretion, Microlife will repair or replace the defective product free of charge.

Opening or altering the device invalidates the guarantee.

The following items are excluded from the guarantee:

- Transport costs and risks of transport.
- Damage caused by incorrect application or non-compliance with the instructions for use.
- Damage caused by leaking batteries.
- Damage caused by accident or misuse.
- Packaging/storage material and instructions for use.
- Regular checks and maintenance (calibration).
- Accessories and wearing parts: Batteries, power adapter (optional).

The cuff is covered by a functional guarantee (bladder tightness) for 2 years.

Should guarantee service be required, please contact the dealer from where the product was purchased, or your local Microlife service. You may contact your local Microlife service through our website: [www.microlife.com/support](http://www.microlife.com/support)

Compensation is limited to the value of the product. The guarantee will be granted if the complete product is returned with the original invoice. Repair or replacement within guarantee does not prolong or renew the guarantee period. The legal claims and rights of consumers are not limited by this guarantee.

## 10. Technical Specifications

---

<b>Operating conditions:</b>	10 - 40 °C / 50 - 104 °F 15 - 90 % relative maximum humidity
<b>Storage conditions:</b>	-20 - +55 °C / -4 - +131 °F 15 - 90 % relative maximum humidity
<b>Weight:</b>	277 g (including batteries)
<b>Dimensions:</b>	131 x 90 x 60.5 mm
<b>Cuff size:</b>	from 17 - 52 cm according to the cuff sizes (see «Selecting the correct cuff»)
<b>Measuring procedure:</b>	oscillometric, corresponding to Korotkoff method: Phase I systolic, Phase V diastolic
<b>Measurement range:</b>	SYS: 60 - 255 mmHg DIA: 40 - 200 mmHg Pulse: 40 - 199 beats per minute
<b>Cuff pressure display range:</b>	0 - 299 mmHg
<b>Resolution:</b>	1 mmHg
<b>Static accuracy:</b>	within $\pm 3$ mmHg
<b>Pulse accuracy:</b>	$\pm 5$ % of the readout value
<b>Voltage source:</b>	<ul style="list-style-type: none"><li>• 4 x 1.5 V alkaline batteries; size AA</li><li>• Mains adapter DC 6V, 600 mA (optional)</li></ul>
<b>Battery lifetime:</b>	approx. 920 measurements (using new batteries)
<b>IP Class:</b>	IP 20
<b>Reference to standards:</b>	IEC 80601-2-30; IEC 60601-1; IEC 60601-1-2 (EMC); IEC 60601-1-11
<b>Expected service life:</b>	Device: 5 years or 10000 measurements, whichever comes first Accessories: 2 years or 5000 measurements, whichever comes first

This device complies with the requirements of the Medical Device Directive 93/42/EEC.

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