
## DESCRIPTION OF THIS PRODUCT GERÄTEBESCHREIBUNG DESCRIPTION DU PRODUIT BESCHRIJVING VAN DIT PRODUCT

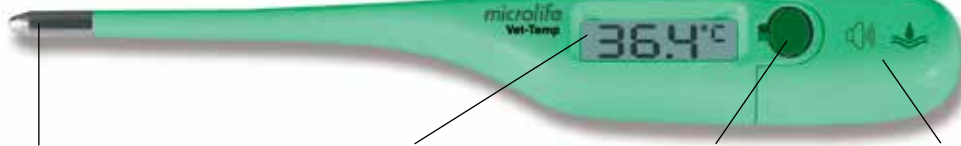Measuring sensor  
Mess-Sensor  
Embout thermosensible  
Meet sensor

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## IMPORTANT SAFETY INSTRUCTIONS

- The thermometer is exclusively for measurement of the animal body temperature.
- The minimum measurement time until the signalling tone [beep] must be maintained without exception!
- Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a minimum distance of 3.3 m from such devices when using this device.
- Protect the thermometer from hitting and striking!
- Avoid ambient temperatures above 60 °C. NEVER boil the thermometer!
- Use only the commercial disinfectants listed in the section «cleaning and disinfecting» to clean the thermometer. The thermometer must be intact when immersed in liquid disinfectant.
- The manufacturer recommends verifying the accuracy by an authorised laboratory every 2 years.

**WARNING:** The measurement result given by this device is not a diagnosis! Do not rely on the measurement result only.

Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

## TURNING ON THE THERMOMETER

To turn on the thermometer, push the ON/OFF button; a short signalling tone [beep] signals «thermometer ON». A display test is performed. All segments should be displayed. Then at an ambient temperature of less than 32 °C, «-» and a flashing «C» appear at the upper right of the display. The thermometer is now ready for use.

## FUNCTION TEST

Proper functioning of the thermometer is tested automatically each time it is turned on. If a malfunction is detected (measurement inaccuracy), this is indicated by «ERR» on the display, and a measurement becomes impossible. In this case, the thermometer must be replaced.

## USING THE THERMOMETER

During a measurement, the current temperature is displayed continuously. The «C» symbol flashes during this period of time. If the signal tone sounds [beep-beep-beep-beep 10 times] and the «C» is no longer flashing, this means that the measured increase in temperature is now less than 0.1 °C in 8 seconds and that the thermometer is ready for reading. To prolong the life of the battery, turn the thermometer off after use by pressing the operating button briefly. Otherwise the thermometer will automatically turn off after about 10 minutes.

## STORAGE OF MEASURED VALUES

If you hold the button down for more than 3 seconds when turning the thermometer on, the maximum temperature stored automatically at the last measurement will be displayed. At the same time, a small «M» for memory will appear at the right of the display. After 3 more seconds, this value will disappear and the thermometer will automatically turn off.

## CLEANING AND DISINFECTING

The thermometer is 100 % waterproof and may be immersed in water or disinfectant solution for cleaning with no problem. Then wipe it off with a moist towel. Please use the following disinfectant for disinfecting, following the manufacturer's instructions:

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**Note:** Any other cleaning solution, dilution or cleaning method can cause a malfunction of the thermometer! The thermometer must never be boiled!

#### MEASURING METHODS

**In the anus (rectal)**  
Carefully insert the tip of the thermometer 2 to 3 cm into the anal aperture.  
Approx. measuring time: 10 seconds!

#### BATTERY REPLACEMENT

When the «▼» symbol (upside-down triangle) appears in the lower right of the display, the battery is exhausted and needs replacing. For the battery replacement please remove the compartment lid of the thermometer. Replace the exhausted battery by positioning the + at the top. Make sure a battery of the same type is on hand. Batteries can be purchased at any electrical goods shops. For battery disposal please refer to «safety instructions».

#### GUARANTEE

We grant you a lifetime guarantee after the date of purchase. Any damage caused by improper handling shall not be covered by the guarantee. Battery and packaging are also excluded from the guarantee. All other damage claims excluded. A guarantee claim must be submitted with the purchase receipt. Please pack your defective thermometer well and send with sufficient postage to:

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